

Media Bottles, Round

Polycarbonate, Non-Sterile



Media Bottles, Round Polycarbonate, Non-Sterile

Our Round Media Bottles are intended for storing and transporting general reagents and cell culture reagents.

Manufactured from durable polycarbonate material, the bottles are supplied without caps. The caps are available separately and are compatible with all three bottle sizes.

The media bottles are chemically compatible with weak acids, alcohols, and aliphatic hydrocarbons. Note that they aren't suitable for exposure to bases, solvents, or other aggressive laboratory chemicals that may cause immediate damage.

The round design minimises contamination with their smoother surface and minimal crevices, allowing for easier cleaning and maintenance. The design also improves fluid circulation and more efficient mixing, in contrast to square bottles where corners can hinder flow. Furthermore, the bottles' structure ensures that pressure is distributed evenly, reducing stress points, and minimising the likelihood of cracks or leaks developing over time.

Product Features

- Media Bottles, Round, Polycarbonate
- Round Design for Reduced Contamination
- Caps Supplied Separately
- Autoclavable, Non-Sterile
- LOT Traceability

Product Specifications

Type	Media Bottles, Round
Material	Polycarbonate (PC)
Colour	Transparent
Bottle Type	Round
Cap Supplied	No
Graduations	Yes



Optional Cap for
Round Media Bottles

Australia
1800 358 101
sales.au@westlab.com
www.westlab.com.au

Canada
1 877 822 1455
sales.ca@westlab.com
www.westlab.com

New Zealand
+64 9 553 3677
sales.nz@westlab.com
www.westlab.co.nz

United States
646 568 5391
sales.us@westlab.com
www.westlab.com

Product Specifications

General Information


LOT Traceability	Yes
Disposable	Yes
Sterility	Non-Sterile
Autoclavable	Yes

Product Information

AUS & NZ	664-371	664-372	664-373
CAN	664-371C	664-372C	664-373C
Capacity	150ml	300ml	600ml
Neck Diameter	37mm	37mm	37mm
Bottom Diameter	50mm	70mm	74mm
Height	117mm	124mm	179mm
Carton Quantity	480	272	168

Product Diagrams


664-371 (150ml)

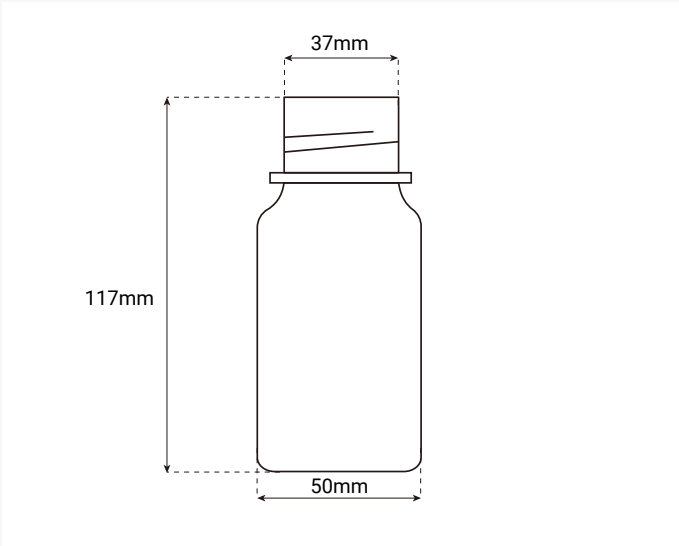


Product Specifications

Product Diagrams

664-372 (300ml)

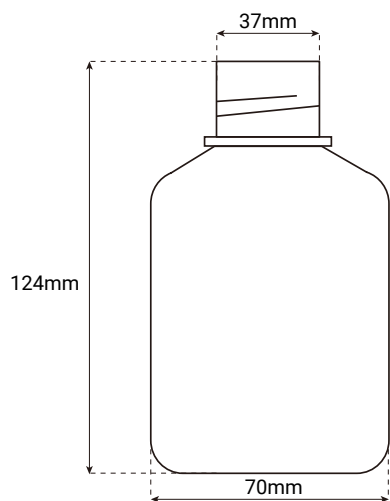

664-373 (600ml)

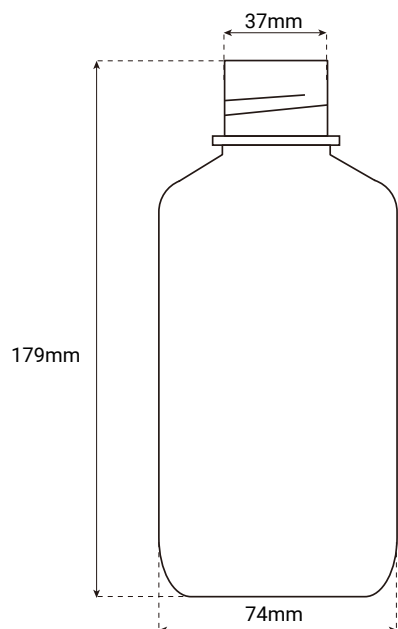

Australia
1800 358 101
sales.au@westlab.com
www.westlab.com.au

Canada
1 877 822 1455
sales.ca@westlab.com
www.westlab.com


New Zealand
+64 9 553 3677
sales.nz@westlab.com
www.westlab.co.nz

United States
646 568 5391
sales.us@westlab.com
www.westlab.com

Optional Accessories

AUS & NZ	664-374
CAN	664-374C
Type	Cap for Round Media Bottles
Material	Polycarbonate (PC) & Silicone Inner Seal
Colour	White
Cap Finish Diameter	37mm
Cap Height	25mm
Sterility	Non-Sterile
Autoclavable	Yes
Carton Quantity	2000

Product Diagram

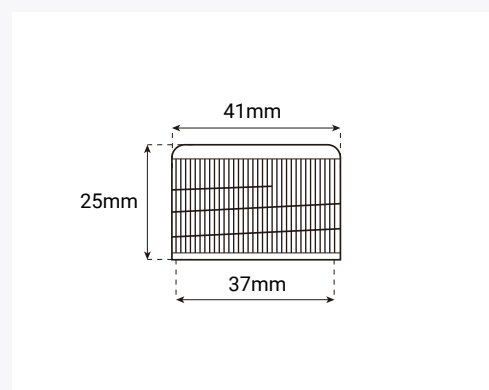

Ordering Information

AUS & NZ	CAN	Description
664-371	664-371C	Media Bottle, PC, 150ml, Round, Non-Sterile, Carton of 480
664-372	664-372C	Media Bottle, PC, 300ml, Round, Non-Sterile, Carton of 272
664-373	664-373C	Media Bottle, PC, 600ml, Round, Non-Sterile, Carton of 168
664-374	664-374C	Cap for Round Media Bottles, PC, Carton of 2000

Publication Date: 08/12/2023

Australia

1800 358 101
sales.au@westlab.com
www.westlab.com.au

Canada

1 877 822 1455
sales.ca@westlab.com
www.westlab.com

New Zealand

+64 9 553 3677
sales.nz@westlab.com
www.westlab.co.nz

United States

646 568 5391
sales.us@westlab.com
www.westlab.com